



www.gmedi.com



EL SERIES


ECG Limb Clamp Electroder



EL-001(Adult)

EL-006(Pediatric)

1 set 4 EA

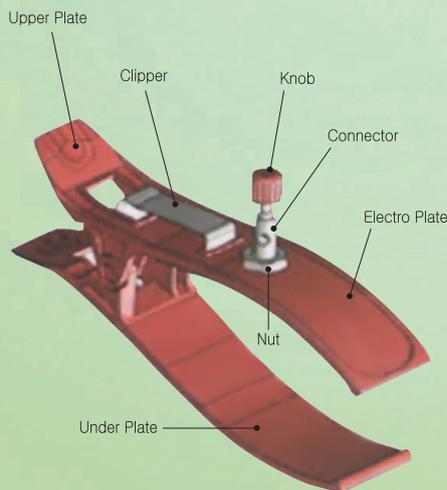

INTRODUCTION

Upper plate Clipper knob Connector Electro plate Under plate nut
EL-001 Electrodes as the accessories of the electrocardiogram, which be applicable in most cardiogram equipments. It could connect with straight Jack, banana Jack and button connector. The Clips are made by special material, for increasing it's elasticity. This product is applicable to all patients (Adult & Pediatric) Ag & AgCl.

ELECTRODES

Remarks : Green, yellow, black and red with universal connector for acquisition of resting E.C.G.

Notes

1. Disinfect electrodes with a 2% solution of a disinfectant commonly used in hospitals.
2. After use, clean electrodes with a soft cloth and a mild detergent. disinfect this products. When use water to clean the products, pls. Attention to the water temp, under 35°C.
3. Do not use during defibrillation!

* The above specifications are subject to change without prior notice.



Manufacture Address:

#401, Chunkang Apartment Shape Factory 402-3, Sindorim-dong, Guro-gu, Seoul, 152-887, Korea
TEL: +82-2-2632-9256 FAX: +82-2-2631-2659 E-mail: gmedir@yahoo.co.kr
<http://www.gmedi.com> <http://eblog.gobizkorea.com/gmedi> <http://kr109859445.trustpass.alibaba.com/>



www.gmedi.com



EL SERIES

ECG Suction Ball Electrode



▲ EL-002
Adult(Electro Plate ø24)

1 set 6 EA

▲ EL-003
Pediatric(Electro Plate ø20)

INTRODUCTION


Knob Suction bowl Electro plate Connector

EL-002 Electrodes as the accessories of the Electrocardiogram. EL-002 Electrodes are fitting for all varieties of ECG Cable & Leadwire. It could connected with straight Jack, banana Jack and snap. Bowl type hardware plated with Ag & AgCl ensures the accuracy and smooth transmittability of Electrocardiograph signal. Which makes electrodes more eligible for the body.

CAUTION

1. Long-time measurements with a strong negative pressure may cause reddening of the skin. Especially when using the electrode on patients with delicate skin, squeeze the suction ball only lightly.
2. When insert into the jack, PLS, tighten the screws, to ensure the jack & are tight interlinked with the Suction Electrodes with no loose.
3. Disinfect the electrodes by wiping them down with a cloth moisten with alcohol (70%).
4. Use available medical cleaning for surgical instruments to clean the electrodes.
5. Do not use metal tools to remove solid particles of dirt as this would destroy the plating.

* The above specifications are subject to change without prior notice.

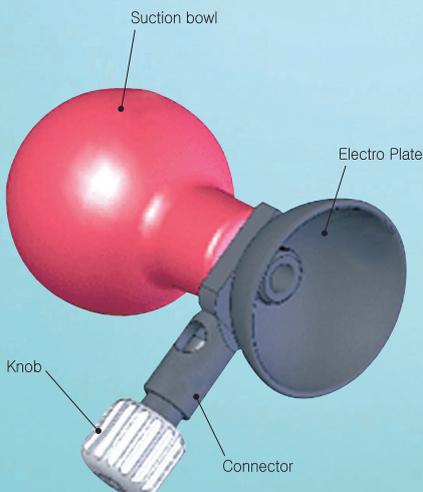

Manufacture Address:

#401, Chunkang Apartment Shape Factory 402-3, Sindorim-dong, Guro-gu, Seoul, 152-887, Korea
TEL: +82-2-2632-9256 FAX: +82-2-2631-2659 E-mail: gmedir@yahoo.co.kr
http://www.gmedi.com http://blog.gobizkorea.com/gmedi http://kr109859445.trustpass.alibaba.com/